





Manual Dialysis Chair **NDC1**



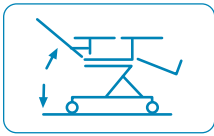

Back Section

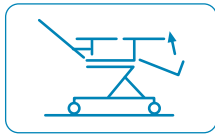

Leg Section

Technical Parameters

• Size	1650mm x 580mm
• Back	770*580mm
• Seat	470*580mm
• Leg	410*580mm
• Sea height	550mm
• Armrest	500x160mm
• Mattress Thickness	60mm

Function

• Back Rest Lifting	20° - 70°
• Leg Rest Lifting	20° - 90°
• Armrest Swing Out	0° - 270°

Technical Configuration

• 3" Silent Caster	1set
• Small Dining Table	1pc

Quality Control

DESIGN
European 3
sections
ergonomics design



STEELS

Baosteel from
Fortune Global 500

PAINTING

11 Process epoxy painting,
ASTM testing anti-bacterial,
paint thickness 0.12mm,
brightness 60°, paint can
resist 50kg impact

0.12_{mm} | 60° | 50_{kg}

WELDING

Panasonic Robotic
ensure
360° full smooth
wedding

360°

ARMREST

Upgraded armrest,
can be right-left,
front-back adjustable.
Donor can control
universal position



Rich Parts



Optional genuine leather or supebric cover, with artisanal marking, comfortable hand feeling.



Adjustable ABS tray table, with antiskid convex edge, tank for cup and cellphone. USB/5V optional. Very convenient for donor.



Full SS IV pole, bearing 15kg





3403 NW 82nd Avenue, Suite 102 E.
Doral, Florida 33122. USA
www.hosmed-inc.com