





Hemodialysis Machine **VANIR**



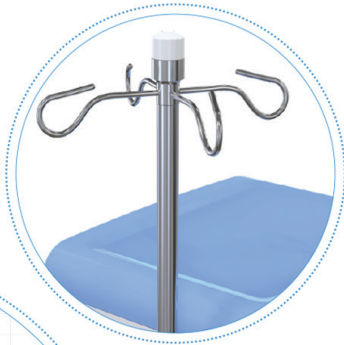

Alarm LED

The IV pole integrates three color alarm LED such as red yellow green, which can be checked from distance with clear version .

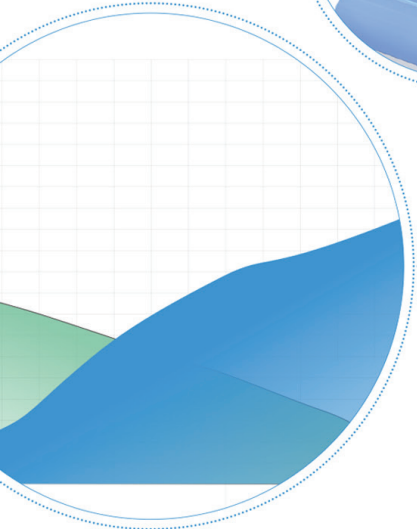


* KTV On-line

The treatment guidance such as Urine clearance, Urea concentration decrease etc can be checked during the treatment, which can provide the treatment prescription.



* BPM

Measure the blood pressure and heart rate of the patients, which can decrease the operation burden of the medical staff.

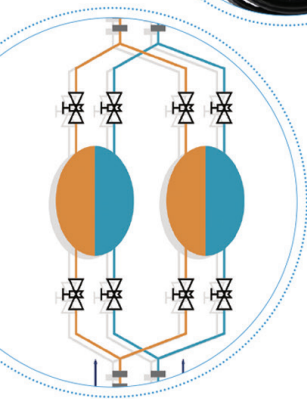

Real-time Leak Monitoring

Balance s system for real-time leak monitoring to ensure the safety of patient.



Back-up Power


Blood circuit runs for more than 40 minutes.






Comprehensive Treatment Mode

- HD
- Sequential HD
- IUF

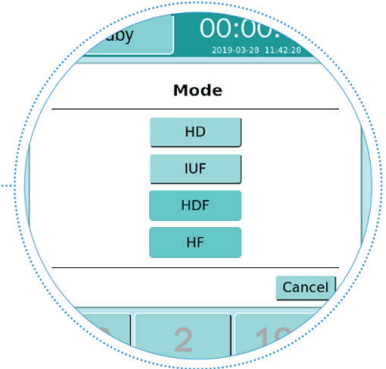

Multi-treatment Profiles, Providing Personalized Treatment

- Sodium profile
- UF rate profile
- Carbonate profile
-
- Heparin flow profile
- Dialysate flow profile
- Dialysate temperature profile

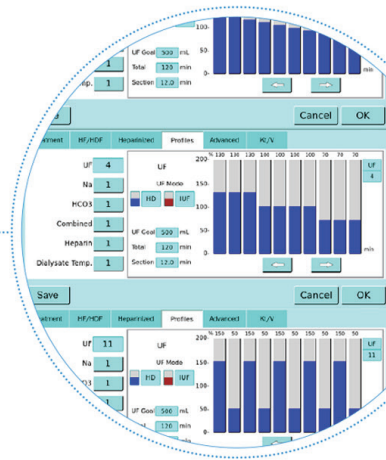

* Bicart Holder

Using Bicart online, avoid microbial contamination of B fluid.




Information System RJ45 Interface

The information management system can be connected through the RJ45 interface, so that realize intelligent and refined management of the entire treatment process.



Technology Empowers Hemodialysis

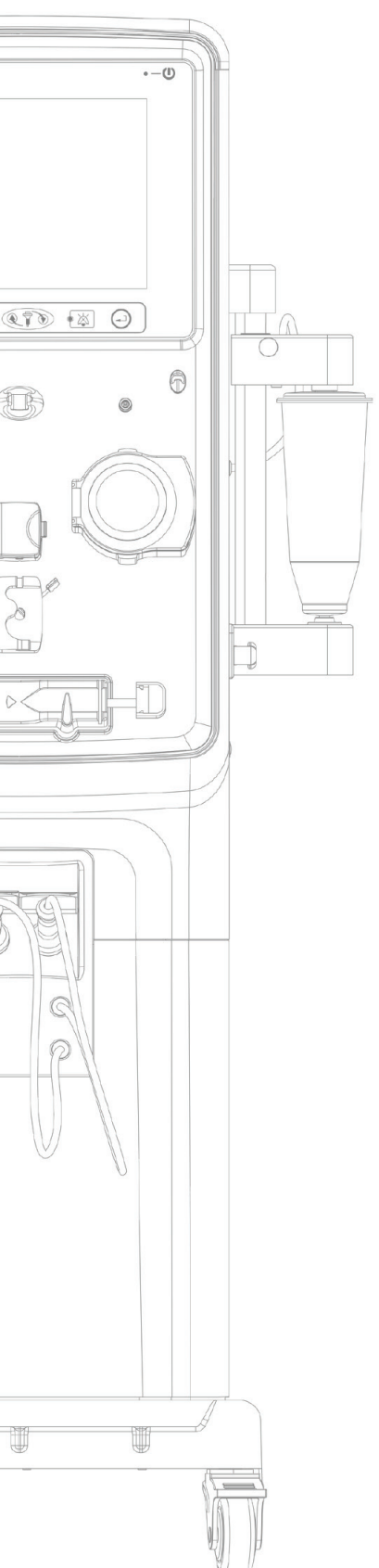

Intelligent Operation

- 12.1 inch LCD+touch screen, multi-language operation interface.
- Adjust the level of venous bubble catcher automatically to prevent blood pollution and reduce dialysis risk.
- Convenient process for sequential dialysis (Dialysis \leftrightarrow IUF) and high-low sodium sequential dialysis.



Safe & Effective Treatment

- The hydraulic, blood circuit and monitor system can be mandatory self-test, so that ensure the safety treatment.
- Ultrasonic and optical detector on venous ensures the safety of blood return.
- Dialysate configuration feedback control system to ensure accurate dialysate concentration.
- Advanced capacity balance feedback control system, ensures UF precision.



Visual Intelligent System

- Preset automatic switch time and disinfection mode, and record all disinfection.
- Dynamic tracking and monitoring functions such as arterial pressure, venous pressure and transmembrane pressure can respond to abnormal situations on time and alarm.



Powerful Compatibility

- Can be operated with Various brands of filters.
- Compatible with various brands of dialysis powder, bloodline and other hemodialysis consumables.



Professional Technical Services

Professional Technical Services: Professional technical services team, which can provide real-time response to hospital requirement and can offer solutions in 12 hours, so that ensure the carefree usage.

Free Regular Training: offer regular training to the medical staffs, which can strengthen sensitivity of medical staffs and emergency response for equipment operation.





3403 NW 82nd Avenue. Suite #102.
Doral, Florida 33122. USA
www.hosmed-inc.com